



BACTRACK® ELEMENT

BREATHALYZER


ESTIMATES BLOOD ALCOHOL CONTENT
IN SECONDS

ENGINEERED FOR MAXIMUM ACCURACY
AND PORTABILITY – IT'S BACTRACK'S
SMALLEST FUEL CELL BREATHALYZER EVER!

Perfect for **Test Yourself At Home™**

- Testing yourself, friends, and family
- Businesses, schools, hospitals, and law enforcement
- Learning how drinking alcohol affects your BAC
- Staying safe and consuming alcohol responsibly

ADVANCED FEATURES Automatic breath flow check and temperature check help ensure proper use.

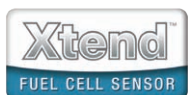

TRUSTED BY PROFESSIONALS Uses police-grade sensor technology for trusted accuracy and consistency.

FAST RESULTS Blow into the mouthpiece and results are displayed in under 10 seconds.

PORTABLE DESIGN Under five inches tall and less than an inch thick – take it anywhere alcohol is consumed.

FOR MORE INFORMATION

© 2011 KHN Solutions LLC. All Right Reserved.



Distributed by: **CLIAwaived.com™**
San Diego, CA 92121
tel 858-481-5031
toll free 888-882-7739
www.cliawaived.com

BACTRACK® ELEMENT

FUEL CELL TECHNOLOGY IN A PORTABLE DESIGN




The BACTRACK Element – Sleek, Portable Design

The new BACTRACK Element quickly and easily tests for the presence of alcohol and provides a Blood Alcohol Content (BAC) estimate. Simply press power, blow through one of the included mouthpieces, and get results in seconds. The Element uses professional Xtend™ fuel cell sensor technology and an internal breath pump to provide trusted results. The Element's small size and affordability make it perfect for anyone – from individuals and businesses to professional groups and even law enforcement.

Test Yourself At Home™

PRODUCT SPECIFICATIONS

Dimensions.....	4.75 x 1.90 x 0.75 inches (12.1 x 4.8 x 1.9 cm)
Weight	3.3 oz (94 g) With mouthpiece and batteries
Sensor Technology.....	Xtend™ fuel cell sensor technology
Detection Range	0.00 – 0.40 %BAC
Mouthpieces	2 included
Test Count	Displays total number of tests performed
Power Supply	Two AAA alkaline batteries, included
Battery Life	Approximately 1000 tests
Warm Up Time.....	10-20 seconds, depending on last use
Response Time.....	10 seconds
Sensor Accuracy	+/- 0.005 %BAC @ 0.1 %BAC
Operating Temperature	32-104 °F (0-40 °C)

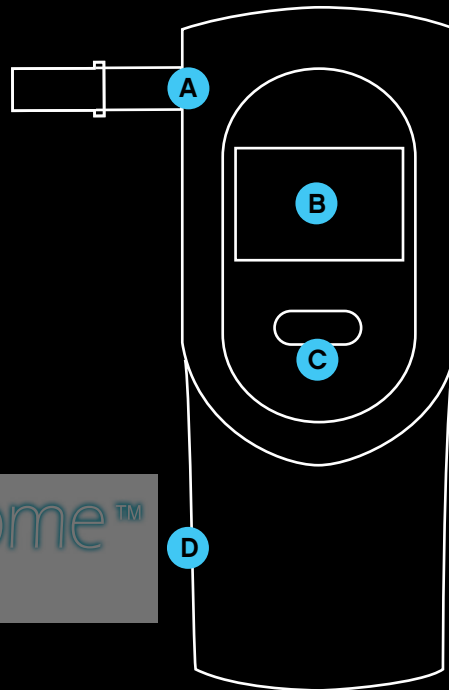

- A** Breath Intake Area – versatile design allows quick use with replaceable, hygienic mouthpieces.
.....
- B** Crisp 4-digit multifunction display – Easy to read test results and alerts.
.....
- C** One Button Operation – Simply turn it on and begin testing in seconds.
.....
- D** Compact size and sleek design - Easily fits into a pocket or purse for on the go lifestyles.

FOR MORE INFORMATION

Distributed by: **CLIAwaived.com™**
San Diego, CA 92121
tel 858-481-5031
toll free 888-882-7739
www.cliawaived.com